



McCallan\*Health. presents

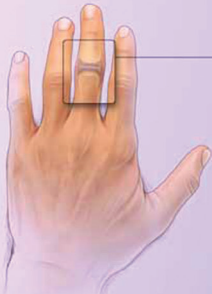

**AMGEN®** *Pfizer*

# Rheumatology Patient Education Tool *Quick Reference Guide*

For Information and illustrative purposes only. Not for diagnostic use.

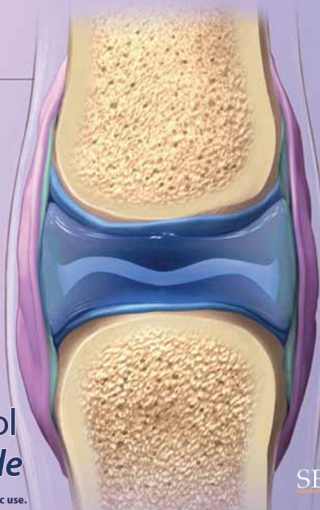

SELECT

## User Icon Key



**Help.** Shows access to Help Desk.



**Settings.** Allows you to choose how to receive updates, etc.



**Information.** Provides privacy and email policies.



**Clear.** Removes the drawing on the image.



**Home.** Takes you back to the image carousel.



**Send.** Allows you to email the drawing or image to the patient.



**Progress Bar.**

Shows how many other images in the carousel.


## Introduction

The Rheumatology Patient Education Tool is provided by Amgen and Pfizer to help you educate your patients about moderate to severe rheumatoid arthritis. It is not intended to diagnose but simply to help illustrate with the use of popular medical and anatomical images. This tool is a web application (app) best viewed on an iPad®.


## How to Use

You will receive access to this app via email. Using your iPad®, click on the URL link provided and your browser will open up to the welcome page.

The first thing you should do is create a link on your home screen to this app. To do this, click on the 'Send' icon on the browser tool bar, and select 'Add to Home Screen'. This will place a Rheumatology icon on your home screen.

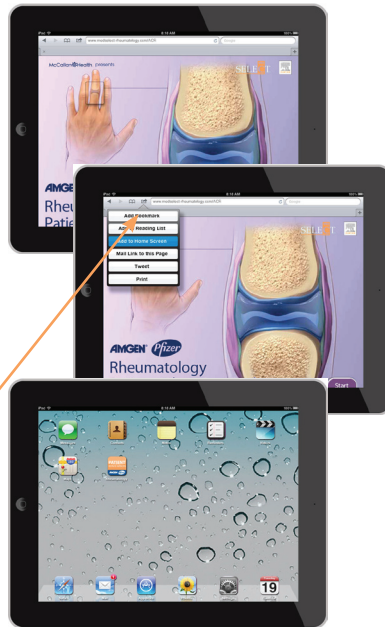

Simply click this icon to launch the Rheumatology app any time you have an internet connection. First time users will need to register and sign in. The screen below shows you this simple step. You should only need to register and sign

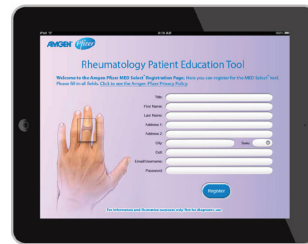

in on the first occasion and will not need to complete this process on subsequent occasions. To ensure that you only need to register the first time you use the app, please set your browser to accept cookies from visited sites. On an iPad® this is found in **Settings > Safari > Privacy**.

Initially you will have been provided a set of images, and additional images will be supplied in the future to qualified physicians. Please take time to familiarize



yourself with these images by scrolling through the image carousel.

When you select an image that you wish to use to illustrate to a patient, simply tap once on the image and that single image will appear.

From here you can pinch in and out to zoom in and out on the image and you can even email this image to your patient (see **Email** section).



Select the drawing icon in the top right hand corner to reveal further options. These options allow you to draw or mark on any given image to enhance the illustration and information you wish to explain to your patient.

The color icons on the left set the line color and the stroke icons on the right set the line thickness.




The Rheumatology Patient Education App can simply be used to illustrate aspects of the disease and treatment to a patient, or you can email them the illustration and what you have marked.

## Email

Please click on the Settings icon to read the email privacy policy which pertains to the Rheumatology Patient Education Tool. In summary, any email is used solely for the purposes of this app and will not be used in any other way, nor will any other communication be sent to the email addresses supplied.

To email an image or drawing to a patient, click on the 'Send' icon. On the following screen (shown below), enter the patient's

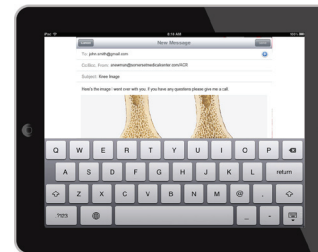

email address and any further information you would like to include. The email will be sent to them via the email address associated with your iPad®.

To email the app development team, please select the *'Help'* icon and click on Help desk.

For any other information concerning Amgen Rheumatology products, please call 1-800-77-AMGEN.

Please visit the information section of the app for details of the privacy policy on email. If you have further questions on email, please contact [help@patienteducation-rheumatology.com](mailto:help@patienteducation-rheumatology.com)

### Notes:

[illegible]



© 2012 Amgen Inc., Thousand Oaks, CA 91320 and Pfizer Inc. All rights reserved. 70510-R1-V1

© 2012 Elsevier Inc. Images from the Netter Collection of Medical Illustrations.

© 2012 McCallan Health LLC. The MED Select logo is a registered trademark and MED Select is a trademark of McCallan Health LLC. All rights reserved.